
Analytical balance KERN ALT-B







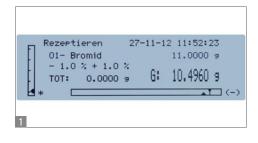

Analytical balance with extensive recipe database

Features

- 11 Convenient recipe-weighing: with the recipe database, in which up to 99 recipes can be stored, each with up to 20 recipe ingredients with name and target value
- Rapid and efficient operation thanks to the graphic display
- · Simple plain text user guidance in the display, following languages available: DE, GB, FR, IT, ES, PT
- Automatic internal adjustment in the case of a change in temperature > 1,5 °C and time-controlled every 20 minutes
- Percentage determination: parts taken out of a container on the weighplate can be displayed as a percentage. Convenient when carrying out drying processes, during which the evaporated moisture or the remaining weight can be displayed as a percentage

- Weighing with tolerance range (checkweighing): Input of an upper/lower limit value. A visual signal assists with portioning, dispensing or grading
- · Vibration-free weighing: When the weighing conditions are unstable, a stable weight is calculated based on an average value
- GLP/ISO record keeping of weight values with date, time and serial number
- Modern all glass draught wind shield for ideal visibility of the weighing object

Technical data

- · Large backlit LCD display, digit height 15 mm
- Dimensions weighing surface Ø 80 mm
- Dimensions housing WxDxH 210x340x330 mm
- Weighing space WxDxH 160x140x205 mm
- Net weight approx. 5,7 kg
- Permissible ambient temperature 5 °C / 35 °C

Accessories

- Protective working cover, KERN ALJ-A01
- Protective dust cover, KERN ABS-A08
- Set for density determination of liquids and solids with density $\leq \geq 1$, for details see page 161. Step-by-step user guidance and calculated density shown immediately on the display, KERN YDB-03
- loniser to neutralise electrostatic charge, see page 162, KERN YBI-01
- Weighing table to absorb vibrations and oscillations, which would otherwise distort the weighing result, see page 162, KERN YPS-03
- · Suitable printers and an extensive accessories range, see page 157 ff.

STANDARD







































Model	Weighing range	Readout	Reproducibility	Linearity	Min. piece weight		Option DAkkS Calibr. Certificate	
KERN	[Max] g	[d] mg	mg	mg	[Counting] mg/piece		DAkkS KERN	
ALT 160-4B	160	0,1	0,1	± 0,3	0,5	0	963-101	
ALT 250-4B	250	0,1	0,1	± 0,3	0,5	0	963-101	

Price reduction

KERN Pictograms



Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven).



Recipe level A: Separate memory for the weight of the tare container and the recipe ingredients (net total).



Suspended weighing: Load support with hook



Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required.



Recipe level B: Internal memory for complete recipes with name and target value of the recipe RECIPE ingredients. User guidance through display.



on the underside of the balance.

Ready for battery operation. The battery type



Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



Recipe level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display. Additional convenient functions, such as barcode and back calculation functions.



Rechargeable battery pack:

is specified for each device.

Rechargeable set.

available.

Battery operation:



Data interface RS-232: To connect the balance to a printer, PC or network.

RS-485 data interface: To connect the balance

tolerance against electromagnetic disturbance.

to a printer, PC or other peripherals. High



Totalising level A: The weights of similar items can be added together and the total can be printed out.



230 V

Power supply: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request.

Mains adapter: 230V/50Hz in standard version

for EU. On request GB, AUS or USA version



RS 485

USB data interface: To connect the balance to a printer, PC or other peripherals.



SUM

Totalising level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display. Additional convenient func-



Strain gauges: Electrical resistor on an elastic deforming body.



Bluetooth data interface: To transfer data from the balance to a printer, PC or other peripherals.



tions, such as barcode and back calculation.



Tuning fork principle: A resonating body is electromagnetically excited, causing it to oscillate.



WLAN data interface: To transfer data from the balance to a printer, PC or other peripherals.



Percentage determination: Determining the deviation in % from the target value (100 %).



Electromagnetic force compensation: Coil inside a permanent magnet. For the most accurate weighings.



Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.



Weighing units: Can be switched to e.g. nonmetric units at the touch of a key. See balance model. Please refer to KERN's website for more



Single cell technology: Advanced version of the force compensation principle with the SC TECH highest level of precision.



Interface for second balance: For direct connection of a second balance.



Weighing with tolerance range: Upper and lower limiting values can be programmed individually for e.g. dosing, sorting and portioning.



Verification possible: The time required for verification is specified in the pictogram.



Network interface: For connecting the scale to an Ethernet network. With KERN products you can use a universal RS-232/LAN converter.



Vibration-free weighing: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value.



DAkkS calibration possible: The time required for DAkkS calibration is shown in days in the pictogram.



GLP/ISO log: The balance displays the weight, date and time, regardless of a printer connec-



Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram. For details see the glossary.



Package shipment: The time required for internal shipping preparations is shown in days in the pictogram.



GLP/ISO log: With weight, date and time. Only with KERN printers, see "Accessories"



ATEX explosion protection: Suitable for use in hazardous industrial environments, in which there is explosion danger. The ATEX marking is specified for each device.



Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram.



Piece counting: Reference quantities selectable. Display can be switched from piece to



Stainless steel:

The balance is protected against corrosion.



Warrantv: The warranty period is shown in the pictogram.

Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2000 kg. In combination with a DAkkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAkkS calibration laboratory today is one of the most modern and best-equipped DAkkS calibration laboratories for balances, test weights and forcemeasurement in Europe.

Thanks to the high level of automation, we can carry out DAkkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- . DAkkS calibration of balances with a maximum load of up to 6 t
- DAkkS calibration of weights in the range of 1 mg 500 kg
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAkkS calibration certificates in the following languages D, GB, F, I, E, NL

Your KERN specialist dealer: